DEUTSCH

ENGLISH

FRANÇAIS

ITALIANO

ESPAÑOL

PORTUGUÊS

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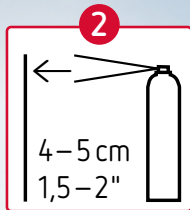
**shake well****spray 1 sec****polish****Lieferform:**

SYMPRO fresh up mint, 75 ml – REF 6700 0300
SYMPRO fresh up lemon, 75 ml – REF 6700 0400

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ENGLISH | INSTRUCTIONS FOR USE

Product description:

SYMPRO fresh up – The finishing touch for cleaned dental acrylic restorations and acrylic surfaces: a satin-matte sheen and a pleasant scent.

Indication:

For the final surface treatment of cleaned dental resin restorations and resin surfaces. For professional use in dental offices and dental laboratories.

Contraindication:

Do not spray directly into the mouth. The product is intended for extraoral use on the resin restoration or surface to be treated. Do not use for cleaning.

Note: This product contains fragrance ingredients that may generally have the potential to cause sensitization.

Special hazard warnings:

Extremely flammable aerosol. The container is under pressure: it may burst when subjected to heat. Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking. Do not spray in the direction of an open flame or other source of ignition. Do not pierce or burn, even after use. Protect against sunlight. Do not expose to temperatures above 50 °C / 122 °F.

Scope of Delivery:

SYMPRO fresh up mint, 75 ml – REF 6700 0300

SYMPRO fresh up lemon, 75 ml – REF 6700 0400

Storage:

- Store at 5 °C to 50 °C / 41 °F to 122 °F.
- Protect against direct sunlight.

Disclaimer:

This product has been developed and is intended exclusively for professional use by trained dental professionals. The instructions for use are to be considered non-binding guidelines. Use and application are the sole responsibility of the user. Renfert assumes no liability for damages that arise from improper handling or use.


Processing:

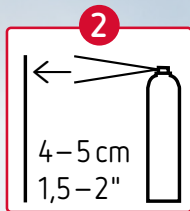- Shake the can well before use.
- Spray the cleaned acrylic surface from a distance of approx. 4–5 cm with a few light spray bursts to ensure even coverage.
- Then distribute evenly with a soft cloth, polish, and wipe away any excess.

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