

The Renfert logo consists of the word "Renfert" in a white, bold, sans-serif font, centered within a red rectangular box. The text is flanked by two white, curved lines that resemble a stylized 'R' or a protective shield.

making work easy

A close-up photograph of a dental technician wearing blue nitrile gloves. The technician is using a dental handpiece with a red, circular, multi-faceted polishing bur to work on a 3D-printed tooth model. The model is held in the technician's left hand, and the handpiece is held in the right hand. The background is a soft, out-of-focus blue and white, suggesting a clinical setting.

# Renfert Polish 3D-printed materials

**NEW** 3D-printing polishing system

# NEW Renfert Polish 3D-printed materials

Polishing specialist for polymer-based printing materials


“Renfert Polish 3D-printed materials” is a polishing paste for 3D-printed resins, designed to smooth layer lines and support marks in a targeted way. It is specifically developed for 3D-printed resins and can be used across different technologies — for example, for filament-printed restorations made of PMMA and PETG as well as for crown and bridge restorations made from resin materials. There is an optional starter kit available, which includes polishing paste and polishers for pre-polishing and high-gloss polishing.

## Advantages

- Specially formulated for polishing 3D-printed temporary restorations.
- Achieve smooth, polished surfaces and a natural shine in no time, thanks to specially selected polishers and a specially formulated polishing paste.
- Ensures hygienic surfaces and comfortable wear by reducing plaque and bacterial buildup.

## Details

- Pre-polishing: Prepolisher 3D side + Prepolisher 3D top
- Polishing: Bison + Renfert Polish 3D-printed materials
- High-gloss polishing: Leather buff + Cotton buff

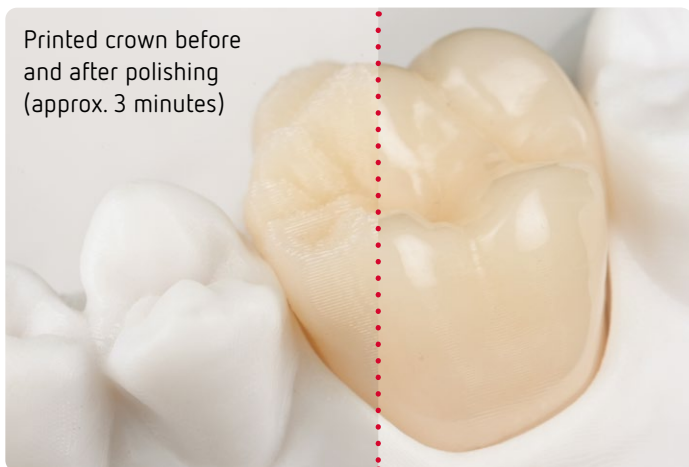

Printed crown before and after polishing (approx. 3 minutes)

## Ordering information

Renfert Polish 3D-printed materials	Polishing paste 13 g	No. 5104000
Prepolisher 3D top	17 mm, 1 piece	No. 880000
Prepolisher 3D side	26 mm, 10 pieces	No. 870000
Bison	14, 18 mm, 25/100 pieces	No. 76xxx00
leather buff	22 mm, 5 pieces	No. 2081000
Cotton buff	22 mm, 12 pieces	No. 2051000
Renfert Polish 3D-printed materials	Starter kit	No. 5104000

